

Choose latex-free catheters.



Universa[®] SILICONE FOLEY CATHETER

Everyday drainage

Choose the Foley catheters that are a part of Cook Medical's full line of products for routine, everyday drainage procedures.

Silicone material

Place a catheter that is comfortable, latex-free, easy to insert, and less prone to plugging and encrustation than latex catheters.¹

Multiple configurations


Select between two- and three-way options, and between various sizes and tip choices for the catheter that fits your preference.

1. Tunney MM, Keane PF, Gorman SP. Assessment of urinary tract biomaterial encrustation using a modified Robbins device continuous flow model. *J Biomed Mater Res.* 1997;38(2):87-93.

Universa[®]

SILICONE FOLEY CATHETER

Universa Silicone Foley Catheters are used to provide continuous bladder irrigation of fluids and drainage of urine from the urinary tract. The catheters can be used for the management of acute retention, obstruction, and intraoperative and postoperative bladder decompression, as well as for urine sampling and monitoring urinary output. They have an indwelling time of up to 30 days.

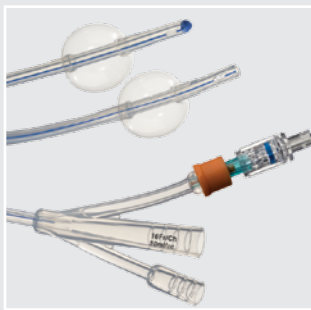

The three-way Foley has lumens for inflation, drainage, and irrigation. Postoperative saline irrigation, for flushing out blood clots in your patients with hematuria, is a typical procedure that requires the use of a three-way catheter.

Order Number	Reference Part Number	Tip	Size Fr	Inflation Volume mL	Working Length cm
Foley Catheter					
Three Way					
G53030	033016-CE	round	16	30	34
G53235	033018-CE	round	18	30	34
G53236	033020-CE	round	20	30	34
G53237	033022-CE	round	22	30	34
G53238	033024-CE	round	24	30	34
G53240	033026-CE	round	26	30	34
G53241	033016-OE	open	16	30	34
G53242	033018-OE	open	18	30	34
G53243	033020-OE	open	20	30	34
G53244	033022-OE	open	22	30	34
G53245	033024-OE	open	24	30	34
G53246	033026-OE	open	26	30	34

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.



The two-way Foley catheter comes in a full range of sizes, from 6-26 Fr. Choose from either a round (closed) tip or an open tip that can be placed over a wire guide. Other options include a larger, 30 mL balloon that is typically used postoperatively or under other special circumstances in which you need a larger balloon size.

Order Number	Reference Part Number	Tip	Size Fr	Inflation Volume mL	Working Length cm
Foley Catheter					
Two Way					
G17111	028506-CE	round	6	1.5	25
G17199	028508-CE	round	8	3	25
G17201	028510-CE	round	10	3	25
G17204	028512-CE	round	12	5	34
G17205	028514-CE	round	14	5	34
G17207	028516-CE	round	16	5	34
G17208	028518-CE	round	18	5	34
G32764	028520-CE	round	20	5	34
G32766	028522-CE	round	22	5	34
G32768	028524-CE	round	24	5	34
G32725	028526-CE	round	26	5	34
G17214	028506-OE-21*	open	6	1.5	25
G17203	028506-OE	open	6	1.5	25
G17200	028508-OE	open	8	3	25
G17202	028510-OE	open	10	3	25
G17198	028512-OE	open	12	5	34
G17206	028514-OE	open	14	5	34
G17345	028516-OE	open	16	5	34
G17512	028518-OE	open	18	5	34
G32765	028520-OE	open	20	5	34
G32767	028522-OE	open	22	5	34
G32769	028524-OE	open	24	5	34
G32771	028526-OE	open	26	5	34
G53230	023016-CE	round	16	30	34
G53229	023018-CE	round	18	30	34
G53231	023020-CE	round	20	30	34
G53232	023022-CE	round	22	30	34
G53233	023024-CE	round	24	30	34
G53234	023026-CE	round	26	30	34

*Includes a .021 inch Slip-Coat® wire guide. Slip-Coat is a registered trademark of Angiotech Biocoatings Corp.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service - Urology

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Poland: +48 223060159, pl.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com
www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
E-mail: cau.custserv@cookmedical.com

